Segmentation of Vestibular Schwannoma from Magnetic Resonance Imaging: An Open Annotated Dataset and Baseline Algorithm (Vestibular-Schwannoma-SEG)

Summary

Redirection Notice

This page will redirect to https://www.cancerimagingarchive.net/collection/vestibular-schwannoma-seg/in about 5 seconds.

This collection contains a labeled dataset of MRI images collected on 242 consecutive patients with vestibular schwannoma (VS) undergoing Gamma Knife stereotactic radiosurgery (GK SRS). The structural images included contrast-enhanced T1-weighted (ceT1) images and high-resolution T2-weighted (hrT2) images. Each imaging dataset is accompanied by the patient's radiation therapy (RT) dataset including the RTDose, RTStructures and RTPlan. Additional ly, registration matrices (.tfm format) and segmentation contour lines (JSON format) are provided and described below. A ll structures were manually segmented in consensus by the treating neurosurgeon and physicist using both the ceT1 and hrT2 images. The value of this collection is to provide clinical image data including fully annotated tumour segmentations to facilitate the development and validation of automated segmentation frameworks. It may also be used for research relating to radiation treatment.

Imaging data from consecutive patients with a single sporadic VS treated with GK SRS between the years of 2012 and 2018 were screened for the study. All adult patients older than 18 years with a single unilateral VS treated with GK SRS were eligible for inclusion in the study, including patients who had previously undergone operative surgical treatment. In total, 248 patients met these initial inclusion criteria including 51 patients who had previously undergone surgery. Patients were only included in the study if their pre-treatment image

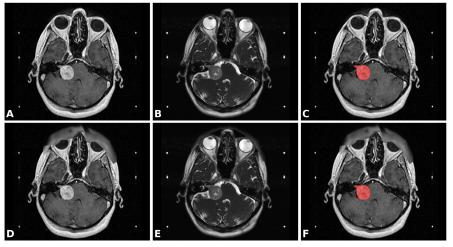


Figure 1. Illustrative example dataset of a patient with a right sided vestibular schwannoma (VS). A) Contrast-enhanced T1-weighted MRI (ceT1); B) High-resolution T2-weighted MRI (hrT2); C) ceT1 MRI with annotated segmentation of VS; D)-F) corresponding images after obscurification of facial features. The six white dots in each image are fiducials of the Leksell Stereotactic System MR Indicator box used for image co-registration.

acquisition dataset was complete; 4 patients were thus excluded because of incomplete datasets and 2 further patients were excluded because the tumour was not completely contained in the image field of view. The images were obtained on a 32-channel Siemens Avanto 1.5T scanner using a Siemens single-channel head coil.

- Contrast-enhanced T1-weighted imaging was performed with an MPRAGE sequence with TR / TE / TI = 1900 ms / 2.97 ms / 1100 ms, in-plane resolution of 0.4 × 0.4 mm, in-plane matrix of 512 × 512, and slice thickness of 1.0–1.5 mm
- High-resolution T2-weighted imaging was performed with a 3D CISS or FIESTA sequence with TR / TE 9.4 ms / 4.23 ms, in-plane resolution of approximately 0.5 × 0.5 mm, an in-plane matrix of 384 × 384 or 448 × 448, and slice thickness of 1.0–1.5 mm.

All manual segmentations were performed using Gamma Knife planning software (Leksell GammaPlan, Elekta, Sweden) that employs an in-plane semiautomated segmentation method. Using this software, each axial slice was manually segmented in turn.

Registration matrices and JSON contours

Registration Matrices: For each subject and each modality there is a text file named inv_T1_LPS_to_T2_LPS.tfm or inv_T2_LPS_to_T1_LPS.tfm. The files specify affine transformation matrices that can be used to co-register the T1 image to the T2 image and vice versa. The file format is a standard format defined by the Insight Toolkit (ITK) library. The matrices are the result of the co-registration of fiducials of the Leksell Stereotactic System MR Indicator box into which the patient's head is fixed during image acquisition. The localization of fiducials and co-registration was performed automatically by the LeksellGammaPlan software.

Contours: For each subject and each modality there is a text file named contours.json. These contour files in the T1 and T2 folder contain the contour points of the segmented structures in JavaScript Object Notation (JSON) format, mapped in the coordinate frames of the T1 image and the T2 image, respectively. There can be small differences between the contour points of the RTSTRUCT and the contour points of the JSON files as explained in the following:

In most cases, the tumour was segmented on the T1 image while the cochlea was typically segmented on the T2 image. This meant that some contour lines (typically for the tumour) were coplanar with the slices of the T1 image while others (typically for the cochlea) were coplanar with T2 slices. After co-registration, the (un-resampled) slices of the T1 and T2 image generally did not coincide; for example, due to different image position and, occasionally, slice thickness.

Therefore, the combined co-registered contour lines were neither jointly coplanar with the T1 nor with the T2 image slices. Upon export of the segmentations in a given target space, the LeksellGammaPlan software interpolates between the original contour lines to create new slice-aligned contour lines in the target image space (T1 or T2). This results in the interpolated slice-aligned contour lines found in the RTSTRUCTs. In contrast, the contours in the JSON files were not interpolated after co-registration, and therefore describe the original (potentially off-target-space-slice) manual segmentation accurately.

Related software

Please see the github respository link which contains a script to organize the downloaded data into a more convenient folder structure and a script to perform co-registration based on the .tfm files and to convert the downloaded DICOM images and JSON contours into NIFTI format. Moreover, the repository contains an algorithm for automatic segmentation of VS with deep learning, adapted to this data set. The applied neural network is based on the 2.5D UNet described in "An artificial intelligence framework for automatic segmentation and volumetry of vestibular schwannomas from contrast-enhanced T1-weighted and high-resolution T2-weighted MRI" (Shapey, et al., 2021) and has been adapted to yield improved segmentation results. Our implementation uses MONAI, a freely available, PyTorch-based framework for deep learning in healthcare imaging (Project MONAI). This new implementation was devised to provide a starting point for researchers interested in automatic segmentation using state-of-the art deep learning frameworks for medical image processing.

Acknowledgements

This work was supported by Wellcome Trust (203145Z/16/Z, 203148/Z/16/Z, WT106882), EPSRC (NS/A000050/1, NS /A000049/1) and MRC (MC_PC_180520) funding. Tom Vercauteren is also supported by a Medtronic/Royal Academy of Engineering Research Chair (RCSRF1819 $\$ \7\34).

Data Access

Data Access

Data Type	Download all or Query/Filter	License
Images and Radiation Therapy Structures (DICOM, 26 GB)	Search Download	CC BY 4.0
Directory names mapping to modality (.csv)	Download	CC BY 4.0

Registration Matrices (.tfm, zip, 257 KB)	Download	CC BY 4.0
Contours (JSON, zip, 16.7 MB)	Download	CC BY 4.0

Click the Versions tab for more info about data releases.

Additional Resources for this Dataset

The NCI Cancer Research Data Commons (CRDC) provides access to additional data and a cloud-based data science infrastructure that connects data sets with analytics tools to allow users to share, integrate, analyze, and visualize cancer research data.

• Imaging Data Commons (IDC) (Imaging Data)

Detailed Description

Detailed Description

Image Statistics	
Modalities	MR, RTSTRUCT, RTDOSE, RTPLAN
Number of Participants	242
Number of Studies	242
Number of Series	1936
Number of Images	48582
Images Size (GB)	26

Citations & Data Usage Policy

Citations & Data Usage Policy

Users must abide by the TCIA Data Usage Policy and Restrictions. Attribution should include references to the following citations:

Data Citation

Shapey, J., Kujawa, A., Dorent, R., Wang, G., Bisdas, S., Dimitriadis, A., Grishchuck, D., Paddick, I., Kitchen, N., Bradford, R., Saeed, S., Ourselin, S., & Vercauteren, T. (2021). Segmentation of Vestibular Schwannoma from Magnetic Resonance Imaging: An Open Annotated Dataset and Baseline Algorithm (version 2) [Data set]. The Cancer Imaging Archive. https://doi.org/10.7937/TCIA.9YTJ-5Q73

(i) Publication Citation

Shapey, J., Kujawa, A., Dorent, R., Wang, G., Dimitriadis, A., Grishchuk, D., Paddick, I., Kitchen, N., Bradford, R., Saeed, S. R., Bisdas, S., Ourselin, S., & Vercauteren, T. (2021). Segmentation of vestibular schwannoma from MRI, an open annotated dataset and baseline algorithm. In Scientific Data (Vol. 8, Issue 1). Springer Science and Business Media LLC. https://doi.org/10.1038/s41597-021-01064-w

(i) TCIA Citation

Clark, K., Vendt, B., Smith, K., Freymann, J., Kirby, J., Koppel, P., Moore, S., Phillips, S., Maffitt, D., Pringle, M., Tarbox, L., & Prior, F. (2013). The Cancer Imaging Archive (TCIA): Maintaining and Operating a Public Information Repository. Journal of Digital Imaging, 26(6), 1045–1057. https://doi.org /10.1007/s10278-013-9622-7

Other Publications Using This Data

TCIA maintains a list of publications which leverage TCIA data. If you have a manuscript you'd like to add please contact TCIA's Helpdesk.

 Shapey, J., Wang, G., Dorent, R., Dimitriadis, A., Li, W., Paddick, I., Kitchen, N., Bisdas, S., Saeed, S. R., Ourselin, S., Bradford, R., & Vercauteren, T. (2021). An artificial intelligence framework for automatic segmentation and volumetry of vestibular schwannomas from contrast-enhanced T1-weighted and highresolution T2-weighted MRI. Journal of Neurosurgery, 134(1), 171–179. https://doi.org/10.3171/2019.9. jns191949

Versions

Version 2 (Current): Updated 2021/03/17

Data Type	Download all or Query/Filter	License
Images and Radiation Therapy Structures (DICOM, 26 GB)	Search Download	CC BY 4.0
Directory names mapping to modality (.csv)	Download	CC BY 4.0
Registration Matrices (.tfm, zip, 257 KB)	Download	CC BY 4.0
Contours (JSON, zip, 16.7 MB)	Download	CC BY 4.0

Registration Matrices (.tfm) and Contours (JSON) were added to the dataset.

Version 1: Updated 2021/02/26

Data Type	Download all or Query/Filter
Images and Radiation Therapy Structures (DICOM, 26 GB)	Search Download
Directory names mapping to modality (.csv)	Download